

# HRC-100 REMOTE CONTROLLER

## Hand Held Laser Remote Controller

The HRC-100 Remote Controller is designed to interface with any Stellar-Pro or Stellar-Pro-L laser system and allows the user to remotely control the function of the laser. The HRC-100 incorporates a number of visual indicators (system power, emission and shutter position indicator LEDs) which allow the user to verify laser status. The HRC-100 also provides manual user control of laser output power as well as shutter position. An 8 foot cable connects the HRC-100 to the laser system via a standard 9 pin D-subminiature male connector.



## Specifications

### Visual Indicators:

System Power Indicator (Red LED)

Emission Indicator (Amber LED)



Shutter Position Indicator (Green LED)

### Laser Controls:

Key Switch Start

Laser Output Power Control

Shutter Open / Close Switch

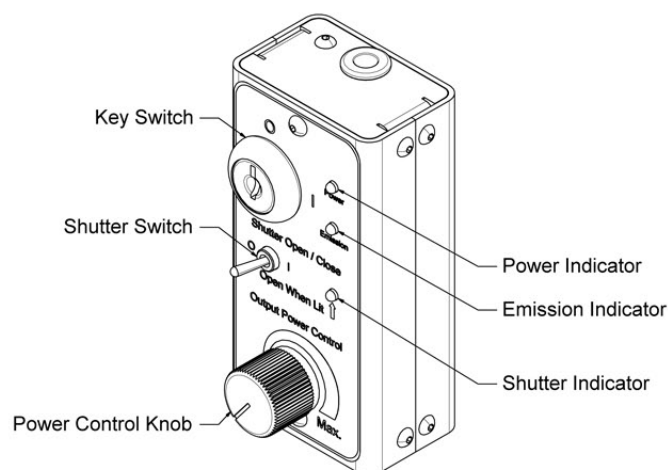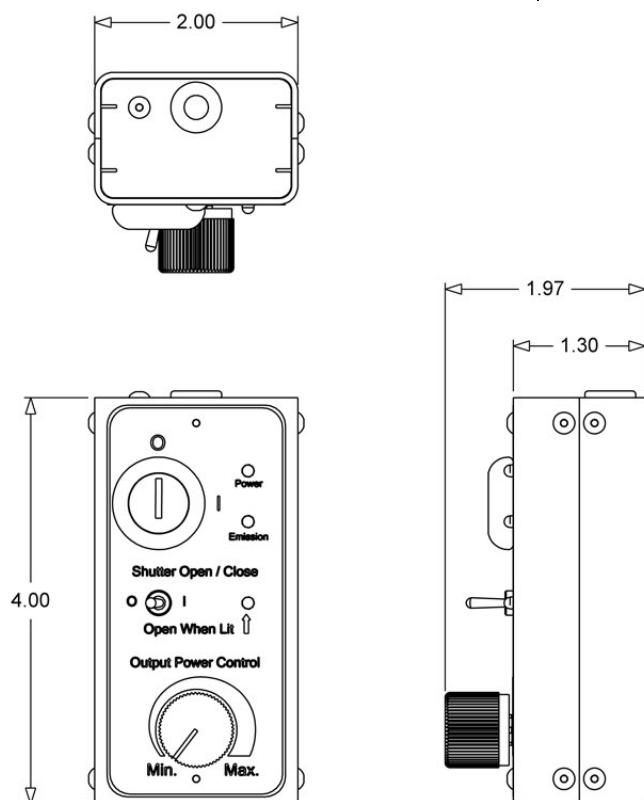

All Dimensions are in Inches

## Modu-Laser, LLC

1258 W. 50 S., Centerville, Utah 84014

Phone: 801-983-2682 Fax: 801-299-2682

e-mail: [info@modu-laser.com](mailto:info@modu-laser.com) web: [www.modu-laser.com](http://www.modu-laser.com)

Specifications are subject to change without notice.

