

RC-1500 LASER CONTROLLER

Argon Laser System Controller

The RC-1500 Laser Controller is designed to interface with Stellar-Pro laser systems (manufactured after June 2012) and provides functional control of the laser system. The RC-1500 incorporates a number of visual indicators (digital laser output power / tube current indicator and system power, emission and regulation indicator LEDs) which allow the user to verify laser status. The RC-1500 also provides manual user control of laser output power as well as shutter control. Connection of the RC-1500 to the Stellar-Pro laser system is achieved via a standard 9 pin D-subminiature (male to male) cable.







Specifications

Visual Indicators:

Digital Laser Power / Laser Current Display
System Power Indicator (Red LED)
Emission Indicator (Amber LED)
Regulation Indicator (Green LED)

Laser Controls:

Key Switch Start
Laser Output Power Control
Shutter "Open" / "Closed" Switch
Operate / Standby Switch

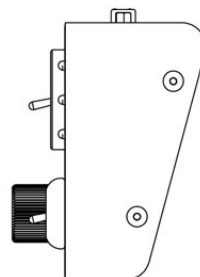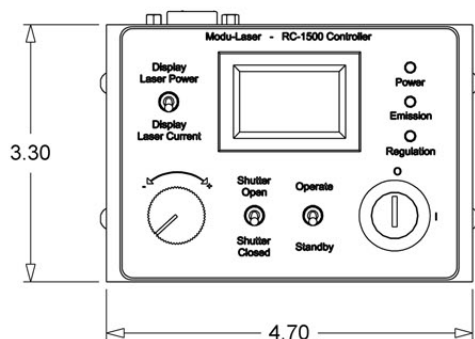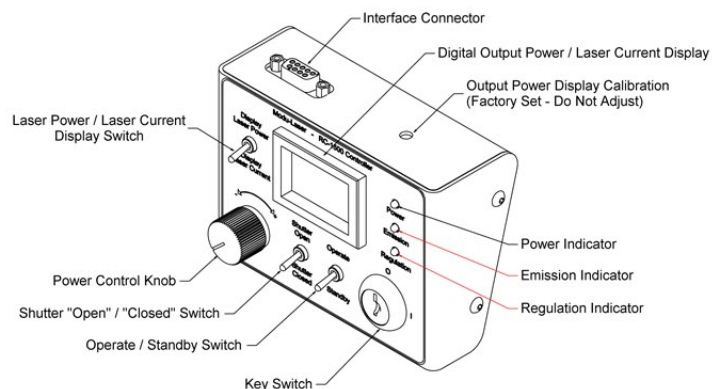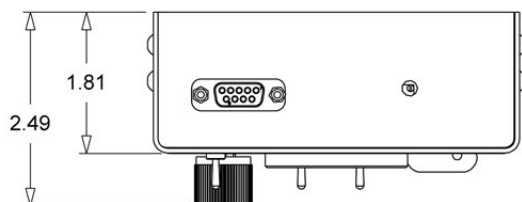

All Dimensions are in Inches

Modu-Laser, LLC

1258 W. 50 S., Centerville, Utah 84014

Phone: 801-983-2682 Fax: 801-299-2682

e-mail: info@modu-laser.com web: www.modu-laser.com

Specifications are subject to change without notice.

